
GETTING SEATED TURNY LOW VEHICLE









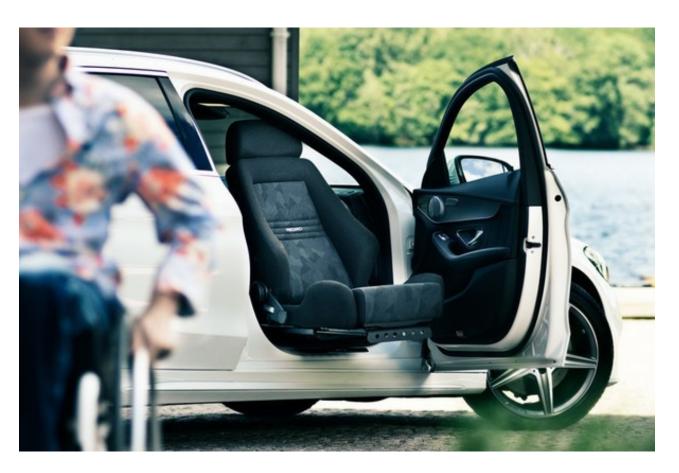

TURNY LOW VEHICLE



Swivels the car seat and extends it almost completely outside the vehicle. Can be installed in many low vehicles previously considered impossible to do.



The **Turny Low Vehicle** is a swivel seat that will allow the user to get seated outside the vehicle, before being comfortably moved to a safe position inside.


Having the car seat outside of the vehicle can be very helpful when getting seated or making the transfer from a wheelchair.

The **Turny Low Vehicle** is the most advanced swivel seat on the market today. Its many unique features makes it possible to adapt car models previously thought impossible to adapt with a swivel seat.



FEATURES

Easy to sit down or transfer from a wheelchair

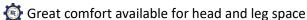



Rotates almost completely outside the vehicle

BENEFITS

Layers can move independently from each other

One layer moves forward, the next layer moves forward, the next layer can rotate and the third can move backwards

Individually customised each movement for you specific need

Can optimise space needed for head, feet and rotation



COMFORT

4 degree inclination to the vehicle floor allowing correct head positioning when exiting the vehicle

Pressure relief

Extra head space and ability to recline the seat without slipping forward





