


GETTING SEATED CARONY 16" (NEW)



GETTING SEATED

CARONY 16" (NEW)

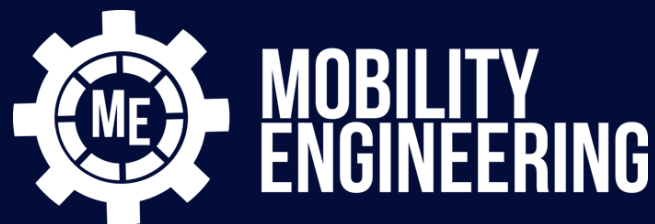

Transfer wheelchair system



The **Carony** is a system that transforms the car seat into a wheelchair and allows for a vehicle transfer that eliminates any manual handling of the occupant. The **Carony** base docks into the swivel seat in the vehicle, either of the Turny variants, and the entire seat easily slides from the car onto the wheelchair base.




The new **Carony 16"** comes with many improvements to help in compatibility with the wider wheelchair market. The occupant experience has vastly changed to help improve comfort and ergonomics as well as additional carer functions such as swing-out handles with brakes.

The seat itself is highly adaptable with individually adjustable armrests, torso support and thigh support. A built in pneumatic bladder allows the front of the cushion to be raised to further improve seating comfort.





This system is the alternative to modifying the vehicle for wheelchair access where the wheelchair occupant does not need a specialised wheelchair.





FEATURES

-  Left or Right applications, integrated push handles with brakes
-  Wheelchair base can be adjusted in height manually
-  Multiple options to cater to occupants seating and postural support

SPECIFICATIONS

-  Base Weight: 15.5 kg + seat
-  Max user weight: 120 kg
-  16" carer assist
-  Crash tested and approved to ISO 7176 and ISO 10542 and CE marked

BENEFITS

-  Functionality of a regular wheelchair with comfort of an automotive seat
-  Suitable for both children and adults providing freedom of independence

