


GETTING SEATED CARONY CLASSIC



GETTING SEATED CARONY CLASSIC

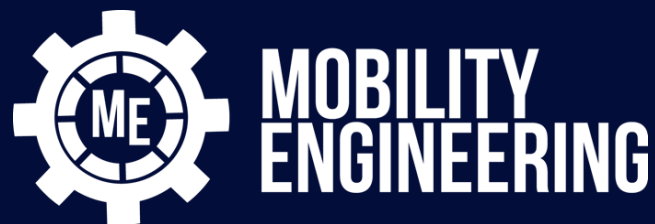

Transfer wheelchair system



The **Carony** is a system that transforms the car seat into a wheelchair and allows for a vehicle transfer that eliminates any manual handling of the occupant. The **Carony** base docks into the swivel seat in the vehicle, either the Turny or Turnout variants, and the entire seat easily slides from the car onto the wheelchair base.




The **Carony Classic** comes in two variants being the 12" and 24" wheelchairs. The 12" **Carony** is a carer assisted wheelchair and the 24" has the option of both carer assist and self-propulsion by the occupant in the wheelchair.

There are 4 seating options for the wheelchair to suit the occupant's needs: **Bev seat, Compact seat, Recaro seat or GS Kids seat**





This system is the alternative to modifying the vehicle for wheelchair access where the wheelchair occupant does not need a specialised wheelchair.



FEATURES



-  Left or Right applications
-  Wheelchair base can be adjusted in height manually
-  Multiple options to cater to occupants seating and postural support

SPECIFICATIONS

-  Weight: 32.5 kg
-  Max user weight: 120 kg
-  12" carer assist or 24" self-propelled model wheelchairs
-  Crash tested and approved to ISO 7176 and ISO 10542 and CE marked



BENEFITS

-  Functionality of a regular wheelchair with comfort of an automotive seat
-  Suitable for both children and adults providing freedom of independence

